

## Exporting ChromaTOF Mass Spectra as .txt File

- 1. Select the compounds you want to export within the desired sample. (To highlight all: highlight first row, hold down SHIFT, then highlight final row.)
- 2. Right click, Export Selected Peaks -> NIST Format...
- Name and save as .MSP file

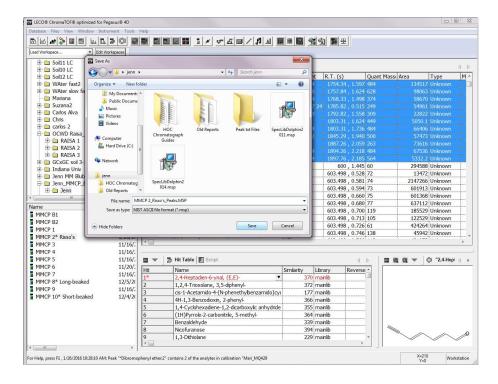

4. Find .MSP file and rename by typing ".txt" as extension

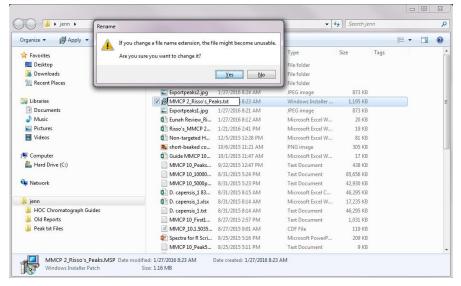

- 5. Select 'Yes'
- 6. Resulting file should appear like so:

MMCP 2_Risso's_Peaks.txt - Notepad	X
File Edit Format View Help	
MRI: Unknown 91.0    Peaks 95   1.5	